


[→ Product Website](#)

**Software for medical professionals**

## Software for calibration and management

RadiCS LE quality control software calibrates EIZO RadiForce monitors and manages the calibration data. Moreover, RadiCS also helps determine and maintain screen data such as inventory numbers and monitor usage time. RadiCS LE features an integrated Backlight Saver function. This allows users to save energy and monitor usage time when the monitor is not in use. RadiCS LE also comes with a Mouse Warp function. This prevents mouse cursors from getting stuck at the edges of screens in multi-screen systems in which the various screens work with different resolutions. RadiCS LE is included with EIZO RadiForce monitors.

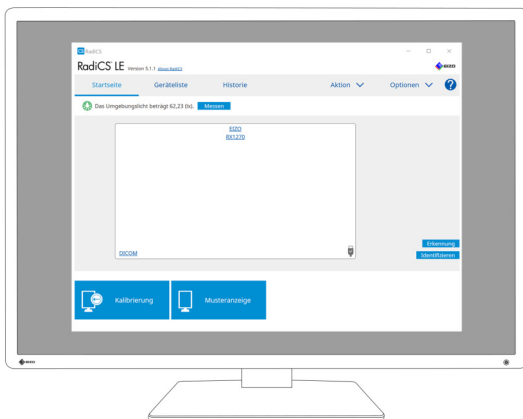

### **RadiCS LE provides you with the following benefits:**

- DICOM calibration of tonal values
- Archives calibration protocols
- Backlight Saver to reduce power consumption and extend service life
- Monitors the sensor inside the monitor that controls brightness and tonal values

## Technical Data

### COMPATIBLE OPERATING SYSTEMS

<b>Windows</b>	Windows 10 / Windows 8.1 / Windows 8 / Windows 7 SP1 / Windows Vista SP2
<b>Mac</b>	OS X Yosemite (10.10) / OS X El Capitan (10.11)

### GENERAL INFORMATION/OPERATION

<b>Display functions</b>	DICOM Part 14 GSDF, CIE, exponentiell (Gammawert), loglinear, linear, benutzerDICOM Part 14 GSDF, CIE, exponentiell (Gammawert), loglinear, linear, benutzerdefinierdefiniert
<b>Interfaces</b>	USB, RS232C, DDC
<b>Languages</b>	Chinese, German, English, Japanese, French

### COMPATIBLE MONITOR

All [RadiForce monitors](#)