


[→ Product Website](#)

Software for medical professionals

Software for Quality control

RadiCS LE quality control software calibrates EIZO RadiForce monitors and manages the calibration data. Moreover, RadiCS also helps determine and maintain screen data such as inventory numbers and monitor usage time. RadiCS LE features an integrated Backlight Saver function. This allows users to save energy and monitor usage time when the monitor is not in use. RadiCS LE also comes with a Mouse Warp function. This prevents mouse cursors from getting stuck at the edges of screens in multi-screen systems in which the various screens work with different resolutions. RadiCS LE is included with EIZO RadiForce monitors.

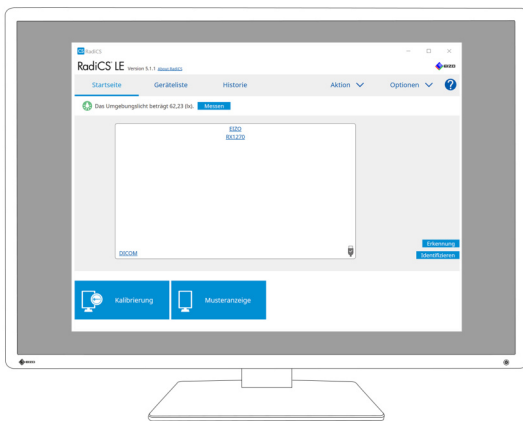

RadiCS LE provides you with the following benefits:

- DICOM calibration of tonal values
- Archives calibration protocols
- Backlight Saver to reduce power consumption and extend service life
- Monitors the sensor inside the monitor that controls brightness and tonal values

Technical Data

COMPATIBLE OPERATING SYSTEMS

Windows	Windows 10 / Windows 8.1 / Windows 8 / Windows 7 SP1 / Windows Vista SP2
Mac	OS X Yosemite (10.10) / OS X El Capitan (10.11)

GENERAL INFORMATION/OPERATION

Display functions	DICOM Part 14 GSDF, CIE, exponentiell (Gammawert), loglinear, linear, benutzerDICOM Part 14 GSDF, CIE, exponentiell (Gammawert), loglinear, linear, benutzerdefinierdefiniert
Interfaces	USB, RS232C, DDC
Languages	Chinese, German, English, Japanese, French

COMPATIBLE MONITOR

All [RadiForce monitors](#)